



Model Number:SDF-001Y
Email address:admin@yanranyl.com
Sales company:Guangzhou Qilai Electronic Technology Co., Ltd
Company address: Room 401, Building G, Jintai Creative Park,
Helong Yiheng Road, Baiyun District, Guangzhou
Company website:www.yanranyl.com

Hengdong Instrument

Electronic massage

Operation Instruction



Credential Information of Hengtong Instrument

Product Name: SDYFU Hengtong Instrument

Product number: SDF-001Y

Attention

1. Please pay attention to personal safety, do not touch plug with wet hands.
2. When using this instrument, please avoid impact and knocking with hard objects to avoid damaging the shell of the instrument.
3. If there is dust on the surface of the instrument, please wipe it with a cloth. Do not use chemicals to clean the instrument.
4. Please pay attention to equal strength when pressing the buttons to ensure long-term normal usage of the instrument.
5. For users with no expressive skills (such as children), adults should be around to watch over them when using the device.
6. Please do not use in rooms with high-temperature (like bathroom) or showering.
7. After use, please make sure that the instrument is turned off, then take out the generator and the massage paste from the main machine, and paste the massage paste back to the translucent plastic sheet. Please do not paste to other things, so as not to lose the viscosity of the massage paste.
8. Please do not optionally press the on-off key to turn off the instrument if the procedure has not completed.

Please read the product instruction and content of Attention before use. Please use the product with the Attentions and other content.

Special Statement

This instruction has been checked and reviewed carefully, but there might exist negligence or mistake of spelling and technology. For such errors, we will make modifications in new instruction without addition notice. The right of interpretation of this instruction belongs to our company, and any other individual or organization shall not revise.

HENG TONG INSTRUMENT

Product display



Front



Back



Side



Bottom

Basic Parameter

Working Voltage: 230V~, 50Hz

Output Voltage: 15V $\overline{\text{---}}$, 0-2.5A

Maximum Power: 60W

Working Temperature: -10° C + 45° C

Product configuration list

Open the outer case of Hengtong Instrument, and check all accessories



1



2



3



4



5



6



7

- 1.Hengtong instrument main engine
- 2.Negative hydrogen ion generator
- 3.Hydrogen molecule generator
- 4.Wrist strap
- 5.Power cord
- 6.Massage patches
- 7.Massage paste power cord
- 8.Instructions



8

It is not suitable to use Heng Tong instrument

Situations of Forbidden to Use Hengtong Instrument

Please turn off the instrument and shut off the power, and contact of after-sale service if following situations happen:

1. Spill water or liquid on the product carelessly, which may cause electric shock.
2. The instrument is damaged and inner parts appear.
3. If you feel abnormal or serious pain when use the product, please stop using immediately and turn to medical advice.
4. If there is any malfunction of machine parts or abnormal situation, please stop using immediately.
5. When current interruption happens in using, please notice that avoid getting hurt from sudden current recovery. When it is lightning and thundering.

Instruction of Operation

Step 1: Please refer to Figure 2. We first plug the power cord into the black jack on the back of the instrument, and then turn on the instrument. Below the jack is the spare fuse (Figure 2/6), the right side (Figure 2/7) is the cooling fan, and the green button (Figure 2/5) is the instrument switch.

Figure 1



1. Negative ion socket
4. Massage socket
7. Heat emission holes

Figure 2



3. Hydrogen molecule socket
- 6.Fuse(include back up fuse inside)

Step 2: Please refer to Figure 1. There are four jacks on the side of the instrument, two generator jacks, which are negative ion socket (Figure 1/1) and hydrogen molecule socket (Figure 1/3), and the wrist strap socket (picture 1/2), massage socket (picture 1/4). Among them, the negative ion generator needs to be used with a wrist strap (please see the product attachment diagram), while the hydrogen molecule generator and massage function can only be used alone. Please select the generator to be used, and then connect it to the corresponding socket on the side of the instrument.

Step 3: Prepare a basin of clean water at about 40°C, and put the connected generator into the water. The height of the water can reach our ankle.

Attention: The wristband should be worn on the inside of the wrist and close to the skin. If skin is very dry, add a little water where the wristband is attached to keep the skin moist.



Figure 3

- 1.Logo
- 5.Time button
- 9.Plus or Minus Display
- 12.Mode Select
- 2.Foot bath Button
- 6 Display Screen
- 10.Time of Foot bath
- 13.Intensity Select

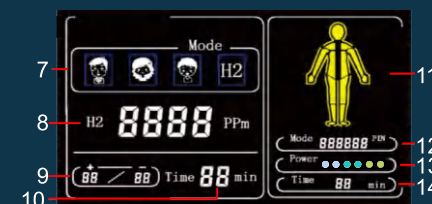


Figure 4

- 3.Set up Button
- 7.Mode Display
- 11.Position of Massage
- 14.Time of Massage
- 4.Massage Button
- 8.PPM Display
- 12.Mode Select
- 13.Intensity Select

Step 4: Foot bath function

The default mode of the instrument is male mode, which lasts for 31 minutes. If you need to choose, please follow the steps below:

1. Mode selection: Tap the foot bath button (Figure 3/2) to choose among four modes: man, woman, child, and H2 (Figure 4/7).
2. Time selection: touch the "+" and "-" buttons (Figure 3/5) to select the foot bath time (Figure 4/10).
3. After confirming the mode and time, press and hold the foot bath button for 3 seconds to start the foot bath.
4. Put your feet in the water and start a foot bath. Please pay attention to the energy value on the display of Hengtong instrument. The normal energy value range is 2500-2800ppb (Figure 4/8). If the displayed value is less than the normal energy value, you can add a small amount of salt to the water.

Step 5: Massage function


The default mode of the instrument is intensity level 1 for 15 minutes. If you need to choose, please follow the steps below:

1. Mode selection: touch and press the button (Figure 3/4), select the function mode from 1-8 (Figure 4/12), (Figure 4-11) is the position of the human body corresponding to the function mode. After selecting the mode, connect the massage stick, and then paste the massage stick on the part we need to massage.
2. Intensity selection: Tap the setting button (Figure 3/3), the intensity display block (Figure 4/13) flashes, tap the "+" and "-" buttons to select the intensity level from 1 to 6.
3. Time selection: Touch the setting button, the time block (Figure 4/14) flashes, touch the "+" and "-" buttons to select the massage time within 15-30 minutes.
4. After confirming the mode, intensity and time, press and hold the massage button for 3 seconds to start working.

Clean and Maintenance Steps of Negative Hydrogen ion Biochemical Generator



1. Hold the generator in one hand and pick up the clamp position (arrow 1) with the finger of the other hand. The cover can be removed.

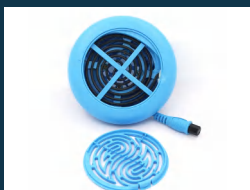

2. The visible part (around the electrode ring) is cleaned with a cleaning brush to remove all the impurities deposited.



3. After cleaning, the positioning point of the cover is aligned with arrow 2, and the position of arrow 1 is pressed with thumb.



4. After cleaning, place in a cool and ventilated place and wait for the next use.

Indications of Hengtong Instrument

1. According to the different constitution of each person, the energy consumption of hydrogen molecular generator is different. When the energy is exhausted, ppm is displayed as 0000.
2. The function of foot soaking and massage can not work simultaneously. To activate one function, the other one will stop.

Corrective Maintenance of Hengtong Instrument

If the product has following malfunctions, please refer to the following relevant methods for solving.

| Malfunction Occurrence | Treatment Methods |
|---|--|
| No display of LED; no electricity of the machine | Check whether the plug is not fix or fuse wire fusing. |
| Hydrogen molecule biochemical reactor does not work | Check whether poor contact of machine or change the hydrogen molecule biochemical reactor. |
| System halted or invalid function of buttons | Turn off the machine and re-start after 10 seconds. |
| The machine has electricity but it does not work | Check whether the metal part of wrist strap touches the skin sufficiently. |