



Hot & Cold Beauty Massager

Instruction manual

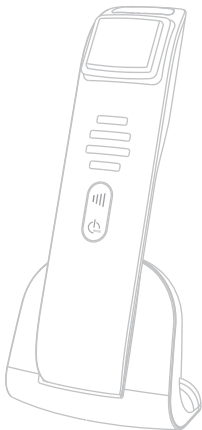

Product structure

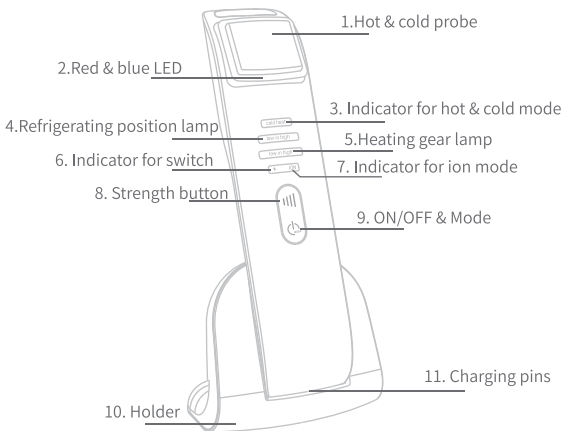

Benefits of hot & cold skincare therapy

Hot & cold skin care therapy can improve skin contractility and enhance skin elasticity. Through 9,000 micro vibrations per minute, it can enhance the vitality of lymphocytes, promote skin metabolism, and improve skin tone. Moreover, it can miniaturize skincare nutrients, improve the penetration into deep skin. Hot therapy gets skin temperature up to 42°C, which is beneficial for the generation of HSP proteins, and dilate pores which is good for deep cleansing. Cold therapy can shrink skin pores, lock moisture in the deep skin, and make skin more delicate and elastic.

The efficacy of cold therapy:

At low temperature 3°C, it can lock moisture deeply, make skin smooth and delicate, reduce intraocular pressure and relieves fatigue for the people with overuse eyes.

At temperature 5°C, it can instantly shrink pores, firm skin, and enhance contours. It's applicable for reducing redness and swelling after aesthetic surgery, also for moisturizing treatment.

At temperature 8°C, it can soothe and calm the skin which suffers from UV sunburn or daily allergies, balance inner skin temperature by cooling down.

The efficacy of hot therapy:

At temperature 42°C, the thermal energy penetrates into deep skin rapidly, activates skin vitality gently, increases skin elasticity, brightens skin tone.

At temperature 45°C, the hot compress improves skin micro-circulation, activates skin, improves dullness, penetrates nutrients into deep skin.

At high temperature 48°C, it's applicable for neck lymphatic drainage massage, the thermal energy helps expanding lymphatic channels, dilates blood vessels and meridians, promotes blood and lymph circulation, discharges heavy metals and toxins.

Sonic vibration:

Sonic vibration massages skin and cells gently, promotes blood and lymph circulation, improves cell permeability & skin metabolism, and makes skin more glossy and elastic.

The efficacy of ions:

Base on the theory of electric charge, which same charges repel but opposite charges attract each other, so the ions+ could lead the dirt out from pores and ions- import nutrient into deep skin.

The efficacy of LED therapy:

Red LED light: penetrate the dermis, activate skin cells and promote the metabolism, rejuvenate skin.

Blue LED light: anti-acne, balance the oil secretion of skin.

Specifications

Power supply: DC3.7(V) 2A	Device net weight: 138g
Battery capacity: 1000mAh	Holder net weight: 58g
Rated power: $\leq 6.5W$	Device size: 182*43*33mm
Charging time: ≤ 3.5 (H)	Holder size: 83*64*77mm

How to use

1. Hold the switch to turn on the device, the default mode is cold mode with Ion- and blue light.
2. Press the strength button to select the strength levels: low-medium-high, loop.
3. Press MODE button to change to heat mode with Ion+ The blue light is off, the red light is on, Hold the strength button to turn on/turn off the vibration.
4. The device will shut down automatically in 5 minutes.
5. The device will beep and shut down automatically while low battery, please charge it.

Charging instruction

Red breathing light on while charging and red light stay on after full charged. The device can't be turned on during charging.

Warranty

1. The warranty period of this product is one year.
2. Please kindly follow the instruction manual to operate the devices, you can get free after-sales services for the failures not damaged by user. Please provide the warranty card together with the faulty product to seller when need maintenance service.
3. Following cases are excluded from the free warranty service even though within warranty period:
 - Failures or damages caused by incorrect operation.
 - Failures or damages caused by disassemble or modify incorrectly by user.
 - Failures or damages caused by the abnormal voltage (adapter with wrong voltage) or natural hazard such as fire, earthquake, flood, lightning stroke.
4. This card will not be issued separately again, please keep it properly.

Hot & Cold Beauty Massager

Warranty Card

Product name	
Purchase location	
Purchase date	
Contact name	
Contact phone	
Delivery address	