



Model Number:SDF-001Y
Email address:admin@yanranyl.com
Sales company:Guangzhou Qilai Electronic Technology Co., Ltd
Company address: Room 401, Building G, Jintai Creative Park,
Helong Yiheng Road, Baiyun District, Guangzhou
Company website:www.yanranyl.com




Hengtong Instrument

Electronic massage

Operation Instruction





Product Name: Hengtong Instrument

Att

Product configuration list

Open the outer case of Hengtong Instrument, and check all accessories



1




2




3



4

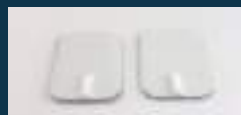


5

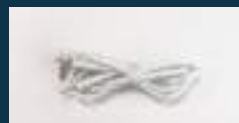

6

- 1.Hengtong Instrument main engine
2. Negative hydrogen ion biochemical generator
3. Wrist Strap
4. Power cord
5. Massage patches
6. Massage paste power cord
7. Instructions

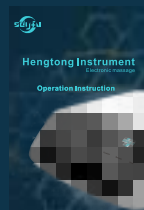

7

It is not suitable to use Heng Tong instrument

Situations of Forbidden to Use Hengtong Instrument

Please turn off the instrument and shut off the power, and contact of after-sale service if following situations happen:

1. Spill water or liquid on the product carelessly, which may cause electric shock.
2. The instrument is damaged and inner parts appear.
3. If you feel abnormal or serious pain when use the product, please stop using immediately and turn to medical advice.
4. If there is any malfunction of machine parts or abnormal situation, please stop using immediately.
5. When current interruption happens in using, please notice that avoid getting hurt from sudden current recovery. When it is lightning and thundering.

Instruction of Operation

In the first step: After we open the case, we should find a place where is near power point and convenient to locate the instrument. Then we can see there is a black socket for power cord, and plug in the power cord. A Green Button which is the on-off key of the instrument.

Picture2



- 1.Socket of generator
- 2.Socket of wrist strap
- 3.Socket of hydrogen ion
- 4.Socket of massage
- 5.Power switch
- 6.Fuse(include back up fuse inside)
- 7.Heat emission holes

Picture1




The second step: As picture 1 shows, there are four sockets at the right side of the instrument, 1.socket for negative ion generator 2. socket for wrist strap 3. socket for hydrogen ion generator 4. socket for massage patches.



Attention: The wrist strap should wear at the inner side of the wrist and touch the skin. If the skin is very dry you can put a little bit of water at the position where skin connected to the wrist strap, to keep skin moist.

The third step: Prepare a basin of water with suitable temprature of foot soaking, normally about 40℃. Put the negative hydrogen ion generator into the water and and position of water surface should 1cm higher then the negative hydrogen ion generator.

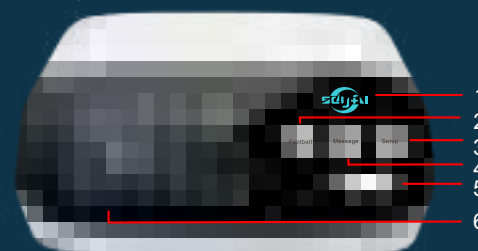

Picture3

- 1.Logo1
- 2.Foot Soaking Button
- 3.Setting Button
- 4.Massage Button
5. Time button
- 6 Display Screen
- 7.Mode Display
- 8.PPM Display
- 9.Plus or Minus Display



Picture4

- 10.Time of Foot Soaking
- 11.Position of Massage
- 12.Mode Select
- 13.Intensity Select
- 14.Time of Massage

The fourth step: To activate function of foot soaking.
The default mode of the instrument is male mode, which lasts 31minutes. You can activate function of foot soaking after 2 seconds of pressing the button of foot soaking. If you need to select, please follow the steps:

- 1.Mode Selecting: to touch the foot soaking button(Male, Female, Kid, H2)four modes for criculation of selectings.
- 2.Time Selecting: touch the time button of “+” and “-”(Picture 3/5)31-60 minutes for selecting.
- 3.After select the mode and time, you can touch the foot soaking button for 2 seconds to activate the function.

The fifth step:
Put feet into water and start to work. Please pay attention to the energy value on the display screen of Hengtong Instrument. The normal energy value ranges from 500-2500. (Picture 4/8) If the displayed number is smaller than the normal energy value, you can add a small number of essence such as cooking salt into the water.

The sixth step: To activate the function of massage
The default mode of the instrument is level 1 of intensity, 15 minutes. You can activate function of massage after 2 seconds of pressing the button of massage. If you need to select, please follow the steps:

- 1. Mode selecting: touch the massage button 1-8 function modes for criculation of selectings.(Picture 4/12)
- 2. Intensity Selecting: touch the setting button (Picture 4/13) flashing, touch the button of “+” and “-”, intensity level from 1 to 6 for selecting.
- 3. Time Selecting: touch the setting(Picture 4/14) flashing, touch the button of “+” and “-”, 15-30 minutes for selecting.

Clean and Maintenance Steps of Negative Hydrogen ion Biochemical Generator



1. Hold the generator in one hand and pick up the clamp position (arrow 1) with the finger of the other hand. The cover can be removed.



2. The visible part (around the electrode ring) is cleaned with a cleaning brush to remove all the impurities deposited.



3. After cleaning, the positioning point of the cover is aligned with arrow 2, and the position of arrow 1 is pressed with thumb.



4. After cleaning, place in a cool and ventilated place and wait for the next use.

Indications of Hengtong Instrument

- 1. According to the different constitution of each person, the energy consumption of hydrogen molecular generator is different. When the energy is exhausted, ppm is displayed as 0000.
- 2.The function of foot soaking and massage can not work simultaneously. To activate one function, the other one will stop.

Corrective Maintenance of Hengtong Instrument

If the product has following malfunctions, please refer to the followin relevant methods for solving.

Malfunction Occurrence	Treatment Methods
No display of LED; no electricity of the machine	Check whether the plug is not fix or fuse wire fusing.
Hydrogen molecule biochemical reactor does not work	Check whether poor contact of machine or change the hydrogen molecule biochemical reactor.
System halted or invalid function of buttons	Turn off the machine and re-start after 10 seconds.
The machine has electricity but it does not work	Check whether the metal part of wrist strap touches the skin sufficiently.